

# Universal Tendon Spacer



## Reconstruction of Flexor or Extensor tendons of the fingers, thumb and wrist

### Characteristics:

- The Universal Tendon Spacer design makes sure that most of the space allowed by the digital canal will be utilized, thus the greatest possible pseudosheath will be obtained.
- The Universal Tendon Spacer can also be used as a physical support to reconstruct the digital canal.
- The Universal Tendon Spacer is made out of High Performance silicone elastomer. It contains a small amount of Bariumsulfate for X-ray viewing.
- The Universal Tendon smooth surface will induce the appearance of a pseudosheath, and in a second step surgery (2 to 6 months later) an active permanent tendon autograft is implanted.

### Dimensions (4 tendon spacers in one)


Length: 49,5 cm

Smaller diameter (B)

equals always half of larger diameter (A)

Diameter A	on a Length of	Diameter B
6 mm – 6 mm	19 cm	3,0 mm
6 mm – 5 mm	8 cm	2,5 mm
5 mm – 4 mm	8 cm	2,0 mm
4 mm – 3 mm	8 cm	1,5 mm
3 mm – 3 mm	6 cm	1,5 mm

*The scale of the drawing is  
1:1 for the width and  
1: ≈2,5 for the length.*

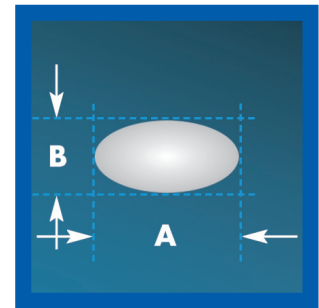

### Manufacturer:

**AREX®**  
*Stay the course*

### AREX®

3, allée du Clos Tonnerre  
BP20  
91125 PALAISEAU Cedex  
FRANCE

### Exclusive Distributor for Germany:



### AMT Aromando Medizintechnik GmbH

Kasernenstraße 21  
D-40213 Düsseldorf  
Telefon: 0211/ 862 999-0  
Telefax: 0211/ 862 999-29  
eMail: info@amt-med.de