www.kimsmed.com

KFDA
식품의약품안전청

Deliver the dream



How is EZ Fan packaged?

Order Code	EZ Fan(REF_01KR)	
	Dimension	Weight
Packaged Device	220 X 330(mm)	36g
Inner Box for 5unit	202 X 277 X 48(mm)	256.11g
Outbox for 100units	515 X 420 X 300(mm)	6.5kg

EZ Fan

Sterile, single use instrument for the delivery of gel implants.



KIMS MED

7 Bonchontaekji-ro, Buk-gu, Gwangju, 61020,
Republic of Korea
T. +82 70 5208 5257 F. +82 70 8611 3022
E. info@kimsmed.com W. www.kimsmed.com

EU REPRESENTATIVE

AR Experts B.V.
Boeingavenue 209, 1119 PD Schiphol-Rijk, The Netherlands

KIMS MED



What is "EZ Fan" ?

EZ Fan is specifically designed to satisfy the needs of the surgeons who want an alternative to the traditional finger insertion method for breast, calf, gluteal and pectoral shaped gel implant procedures.

What is the indication for "EZ Fan" ?

EZ Fan is indicated to deliver gel implants into the pocket with minimal or no touch of the implant, surrounding tissue and the incision.

What are the advantages from "EZ Fan" ?

- Employ a No-Touch / Minimal-Touch technique believed by respondents to reduce infection rates
- Shorten incision lengths
- Reduce operating room time
- Reduce force applied to the implant shell, which may lessen damage occurrence
- Increase placement options
- Can accommodate all types of gel implants
- Less friction, faster propulsion of implant
- Superior visualization of implant

How to use "EZ Fan" ?

HYDRATE


LOAD

TRIM

PROPEL

HYDRATE

Fully immerse the funnel in a basin of sterile solution to ensure that the inner surface is slippery. When not in use, HydroCone should be laid flat on a sterile table so that the inner surface remains moist. The inner surface should be fully hydrated before using.

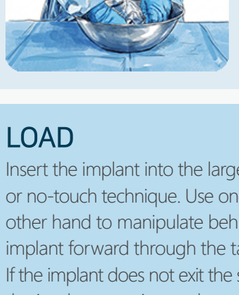

LOAD

Insert the implant into the large proximal end of the hydrated funnel using a minimal or no-touch technique. Use one hand to support the distal end of the funnel, and the other hand to manipulate behind the implant. Use this hand movement to push the implant forward through the tapered funnel and eject the implant into a sterile field. If the implant does not exit the small end of the funnel using gentle force, then remove the implant, re-trim, and repeat the test, gradually trimming larger opening sizes until the desired result is achieved. Load the implant into the funnel ready for delivery into the tissue pocket. Using one hand, apply slight pressure behind the implant to remove any slack, verify positioning of the implant orientation within the funnel and advance the implant to approximately 1cm from the distal end.




TRIM



It is recommended to use sterile straight edge blade scissors to trim the distal end. The HydroCone opening should be large enough to allow the implant to gently pass through without being too large, thus inadvertently enabling the implant to pass through.

Information to trim the EZ Fan funnel	
Smooth Round	Textured Round
Moderate&Moderate Plus : as shown	Moderate&Moderate Plus : +1.5cm larger
High Profile: +0.5cm larger	High Profile: +2cm larger
Ultra High Profile: +1cm larger	Ultra High Profile: +2cm larger
Anatomical / Shaped Implant +2cm larger	

PROPEL

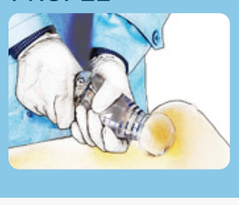

Retract the incision site fully and continue to apply slight pressure behind the implant. Insert the distal end of the funnel approximately 1 cm through the incision. Use one hand to support the funnel and use the opposite hand to employ a pushing/squeezing technique behind the implant to gently propel the implant forward into the centre of the surgical pocket.